

Clinispin horizon 6, 12 & 24 Flex

Programmable Routine Centrifuges

Satisfy the most diverse processing requirements with Clinispin horizon Flex Routine Centrifuges. Customise settings and check every detail on the digital display or pre-program up to 10 cycles.



- Fully Programmable
- Rapid Horizontal Separation
- 2 Stage Horizontal Rotor
- Clear Shatterproof Lid for Easy Monitoring
- Lid Safety Interlock System
- Brushless DC Motor
- Cool Operation

Clinispin horizon 6 Flex

Ideal for small routine labs spinning a variety of samples or seeking a centrifuge with a fully digital display

Clinispin horizon 12 Flex

Ideal for medium sized labs with a variety of sample preparations or where a digital display is desirable

Clinispin horizon 24 Flex

Ideal for large labs with many varied specimen preparations or needing the second cycle countdown

Suggested Applications

- Doctor's Offices
- Small to Large Laboratories
- University Departments
- Multiple Tube-Size Requirements
- Contract Research/
Pharmaceutical Organisations

Visit www.woodleylabdiagnostics.com for full product range

Clinispin horizon 6, 12 & 24 Flex Programmable Routine Centrifuges

Satisfy the most diverse processing requirements with Clinispin horizon Flex Routine Centrifuges. Customise settings and check every detail on the digital display or pre-program up to 10 cycles.

Lid lights indicate cycle status:
ready, running, complete

Display shows remaining time
Verify exact cycle parameters



Use default presets* for chemistry, coagulation and urine or program your validated cycles


*presets based on CLSI guidelines and tube manufacturer IFU

10 programmable settings for all your applications

Features

- **Fully Programmable** - save & recall 10 easily customisable cycle settings or choose from 3 presets
- **Lid Lights** - flash at the end of each cycle & continue flashing until the lid is opened
- **Automatic Lid Lock** - will not run with lid open & lid cannot be opened while samples are spinning
- **Cool Operation** - unique design of the convection based cooling system protects the samples even under the heaviest workloads

- **Clear Shatterproof Lid** - for easy monitoring
- **Brushless DC Motor** - no routine maintenance
- **Reinforced Guard Bowl** - maximum safety & durability
- **2 Stage Horizontal Rotor** - patented swing-out design produces ideal horizontally separated samples

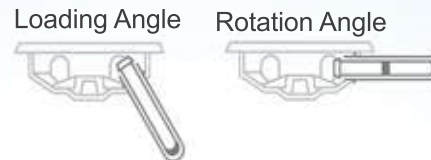

Rotor and Tube Holders Included

Specifications

	HORIZON 6 Flex	HORIZON 12 Flex	HORIZON 24 Flex
Tube Capacity	6 tubes – 3-10 mL	12 tubes – 3-10 mL	24 tubes – 3-10 mL 12 tubes – 3-15 mL 6 tubes – 50 mL
Dimensions (H x W x D)	36 x 30 x 23 cm (14" x 12" x 9")	33 x 38 x 23 cm (13" x 15" x 9")	38 x 43 x 23 cm (15" x 17" x 9")
Weight	5.4kg (12 lbs)	15kg (34lbs)	17kg (39 lbs)
Sound Level	64 dB A	56 dB A	59 dB A
Environmental Range	16-32°C	16-32°C	16-32°C
Voltage	95-253 VAC	95-253 VAC	95-253 VAC
Frequency	50/60 Hz	50/60 Hz	50/60 Hz
Power Requirement	220 Watts	415 Watts	415 Watts
Centrifuge Motor	½ H.P. Brushless	½ H.P. Brushless	½ H.P. Brushless
Max G-force	2000 xg	2000 xg	2000 xg
Max Speed	3800 RPM	3700 RPM	3400 RPM
Cycle Time	1-30 minutes (+/- 2%)	0.5 to 99 minutes (+/- 2%)	0.5 to 99 minutes (+/- 2%)

Visit www.woodleylabdiagnostics.com for full product range