


DiaSys Blood Urea Test Strips

Test strips for the determination of urea
in whole blood, serum and plasma

A quick and easy, low cost test for the quantitative determination of urea in-practice. This enzymatic assay utilises urease to split urea into carbon dioxide and ammonia. Ammonia colours the test strip blue and migrates up the strip scale. The length of the blue zone indicates the urea concentration providing an immediate assessment of the patient's renal function.

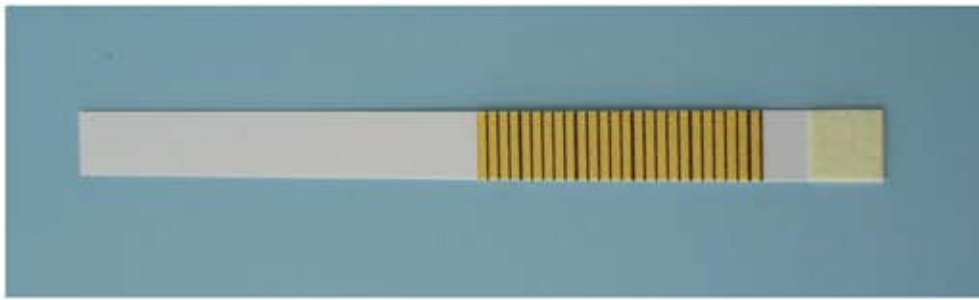

Strip Storage

Room temperature (15-25°C)

Sample

10µl capillary, whole blood, serum or plasma

Pack size

60 Strips

Test Time

From 5 minutes

Measuring Range

5 - 64 mmol/L

14 - 180mg/dL



*Quick and easy
assessment of renal
function.*



Tel: +44 (0)1204 669033

Fax: +44(0)1204 669034

Email: sales@woodleyequipment.com

Web: www.woodleyequipment.com



WOODLEY
EQUIPMENT COMPANY LTD.

OLD STATION PARK BUILDINGS, ST. JOHNS STREET, HORWICH, BOLTON, LANCASHIRE, BL6 7NY, UK