

InSight[®] Urinalysis

Available 24 hours a day



Results in 15 minutes



Results within 24 hours

IN-PRACTICE URINALYSIS.

CONFIRMATION OF UTI IN 15 MINUTES.

ANTIBIOTIC SENSITIVITY WITHIN 24 HOURS.

NO NEED FOR SPECIALIST KNOWLEDGE.



"PASSIONATE ABOUT LABORATORY DIAGNOSTICS"

WOODLEY
EQUIPMENT COMPANY LTD.

InSight[®] Urinalysis


A 2-stage process available 24 hours a day, with confirmation of UTI in 15 minutes.

Perform in-house without the need for specialist knowledge

1. RapidBac[™] Vet 2.0

TEST FOR BACTERIURIA

- * Negative - no bacteriuria indicated
- * Positive - proceed to Flexicult Vet if antibiotic sensitivity required

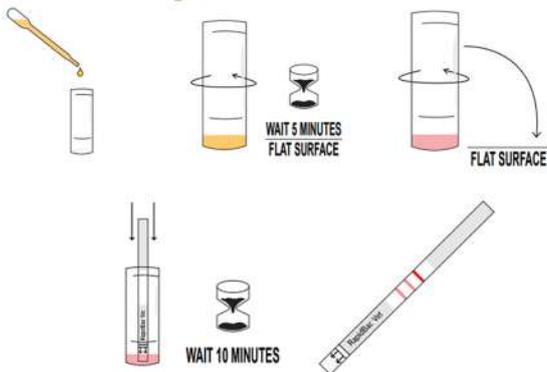

Results in 15 minutes

RapidBac Vet 2.0 is a lateral flow immunoassay for the detection of Gram-positive and Gram-negative bacteria in Dog and Cat urine.

- Bufferless formation
- Test immediately in-practice
- Results in 15 minutes
- Easy to use
- Room temperature storage
- 98.5% accuracy
- No special apparatus required

Proceed with Flexicult Vet if RapidBac is positive & antibiotic sensitivity is required

2. Flexicult[®] Vet

Chromogenic In-Practice Urine Culture & Sensitivity Test



Flexicult Vet is a 24 hour chromogenic culture test for diagnosing UTI's in Dogs and Cats. It aids the identification and quantification of bacteria and helps select the most appropriate antibiotic for effective patient management.

Results within 24 hours

Flexicult Vet contains the following antibiotics:

- Ampicillin
- Amoxicillin/clavulanate
- Oxacillin
- Enrofloxacin
- Trimethoprim/sulfamethoxazole

Agar plate contains a chromogenic substance

Causes bacterial colonies on the agar to appear with different colours depending on the bacterium species



Examples of Flexicult Vet inoculated with E. coli (left) and E. faecalis (right)



Incubator required - available from Woodley

+44 (0) 1204 669033 (Option 1)

sales@woodleyequipment.com

www.woodleyequipment.com

Old Station Park Buildings
St. John Street
Horwich
Bolton
Lancashire
BL6 7NY
UK

"PASSIONATE ABOUT LABORATORY DIAGNOSTICS"

WOODLEY
EQUIPMENT COMPANY LTD.