

InSight V5[®]R

Veterinary Haematology Analyser



**5 PART DIFFERENTIAL
WITH RETICULOCYTES.**

"PASSIONATE ABOUT LABORATORY DIAGNOSTICS"

WOODLEY
EQUIPMENT COMPANY LTD.

InSight V5[®]R

A fully automated multi-species,
5 Part Differential Haematology
Analyser with Reticulocytes

FEATURES

- Quick and easy to use
- Small sample volume - 20µl whole blood
- 5-part differential, 25 parameters
- Reticulocytes can be selected on every sample or selectively on each sample
- Species specific reference ranges
- Independent channel for Basophil measurement
- Powerful capability to flag abnormal cells
- Fully automated
- Touch screen
- Low maintenance
- Results transmitted via InSight p-LINK to Practice Management Systems

SPECIFICATIONS

Principles

- Flow cytometry & Tri-anglescatter for 5-part differential analysis and WBC counting
- Impedance method for WBC, RBC and PLT counting
- Cyanide free reagent and colourimetric method for haemoglobin concentration testing
- Laser scattering method and staining method for Reticulocytes

Parameters

WBC	LYM#	HGB	HCT	Retic#
NEU%	EOS%	MCV	PLT	Retic%
NEU#	EOS#	MCH	MPV	
MON%	BASO%	MCHC	PDW	
MON#	BASO#	RDW-CV	PCT	
LYM%	RBC	RDW-SD		

Species:

- Group I: Dog, Cat
- Group II: Rabbit, Cow, Horse, Sheep, Pig
- Group III: Mouse, Rat, Guinea Pig, Monkey, Camel

User Interface


- Easy to use, menu driven user interface
- Screen Type - LED backlit, colour graphics LCD

Connection Interfaces

- One LAN interface
- Printer Options - Built-in thermal printer or connect to external inkjet and laser printers
- 4 x USB interfaces



Do you need retics checking?
Save on reagents for each sample

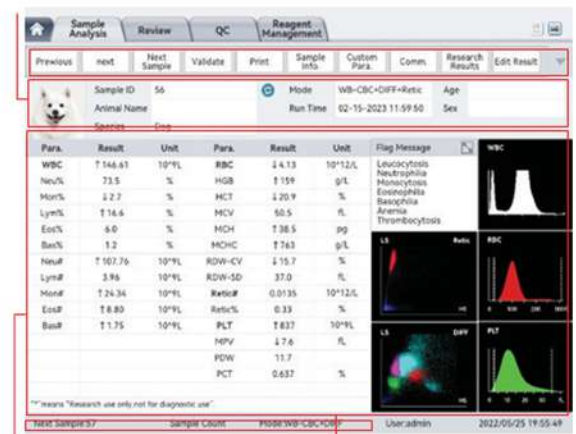

INTEGRATES WITH PRACTICE MANAGEMENT SYSTEMS

Graphic Parameters

- Five 2D scattergrams (including 1 for Baso, 1 for Retics)
- Three histograms (including RBC, PLT histograms)

Software Languages

English, Chinese, Spanish, French, Russian, Portuguese

Dimensions

- Width 36cm (14.3")
- Depth 50cm (19.7")
- Height 49cm (19.3")
- Weight 32kg (70.5lbs)
- Power Supply - 100-240V AC

+44 (0) 1204 669033 (Option 1)

sales@woodleyequipment.com

www.woodleyequipment.com

Old Station Park Buildings
St. John Street
Horwich
Bolton
Lancashire
BL6 7NY
UK

"PASSIONATE ABOUT
LABORATORY DIAGNOSTICS"

WOODLEY
EQUIPMENT COMPANY LTD.