


www.woodleyequipment.com

A stand alone Lab Management System to integrate with your Practice Management System

InSight™

Lab Manager

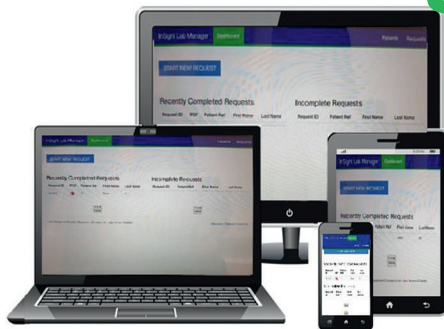

Results available on your smartphone or tablet - access Lab Results anytime, anywhere

Compatible with

BIONOTE Vcheck V200

Achieve true analysis with fluorescence technology

- Canine Cortisol, Canine & Feline T4

- Pancreatitis Markers cPL, fPL

- Real-Time Inflammation Markers cCRP, fSAA

see page 4 for more information

Wet Chemistry, Endocrinology and Blood Gas Options



- InSight Lab Manager integrates your Lab Equipment with the host Practice Management System
- Alternatively InSight Lab Manager can be used as a stand alone Lab Management System
- Results are displayed alongside a Reference Range Bar Chart for easy identification of abnormal results

Contact Woodley or your Local Area Account Manager to discuss

Tel: +44 (0) 1204 669033 (Option 1)

Fax: +44 (0) 1204 669034

Email: sales@woodleyequipment.com

Web: www.woodleyequipment.com

OLD STATION PARK BUILDINGS, ST. JOHN STREET, HORWICH, BOLTON, LANCASHIRE, BL6 7NY, UK


WOODLEY
EQUIPMENT COMPANY LTD.

A stand alone Lab Management System to integrate with your Practice Management System

InSight™ Lab Manager

Results available on your smartphone or tablet - access Lab Results anytime, anywhere

- InSight Lab Manager operates on a secure isolated network within the Practice Laboratory
- InSight Lab Manager is designed from the bottom up to be smartphone and tablet friendly
- Combining the InSight Lab Manager application with a 3G connection, users can access the application functions anytime, anywhere
- The results can be automatically transmitted to the host Practice Management System or printed directly from the InSight Lab Manager system on any networked printer
- Results are displayed alongside a Reference Range Bar chart for easy identification of abnormal results

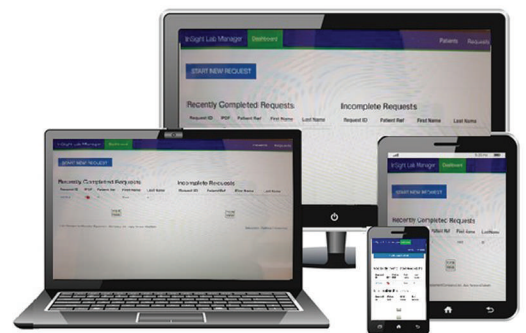

How does it work?

- InSight Lab Manager connects to each analyser using a secure isolated network
- InSight Lab Manager connects to the Practice Management System via LAN or WiFi or can be accessed via 3G if required
- The user starts a lab 'request' on InSight Lab Manager which will then generate a unique request ID
- The samples are processed on the analyser(s)
- InSight Lab Manager then receives the results and matches them to the request
- A **Lab Report** is then autonomously compiled and transferred to the Practice Management System

WOODLEY EQUIPMENT COMPANY LTD						
Woodley Equipment Company Ltd Old Station Park Buildings, St. John Street, Horwich, Bolton, BL6 7NY						
Laboratory Report						
Patient Detail						
Name	Ref	Sex	Age	Species	Breed	
Bob Shaw	18093	male	5	canine	Boston terrier	
Results Detail						
Test	Result	Unit	Range	Flag		
ALB	31	g/L	20-45			
TPP	73	g/L	50-80			
GLU	5.3	mmol/L	3.3-6.5			
ALP	120	U/L	0-212			
ALT	53	U/L	0-48			
BUN	7.1	mmol/L	2.1-6.3			
CRISA	106	mmol/L	35-145			
WGLDB	42	g/L	25-45			
BMG	0.7		0.0-0.0			
BVC	47.0		0.0-0.0			

All test results from either lab analysers or reference laboratories can be captured and assimilated onto a single pdf lab report

Tel: +44 (0) 1204 669033 (Option 1)
Fax: +44 (0) 1204 669034
Email: sales@woodleyequipment.com
Web: www.woodleyequipment.com

WOODLEY
EQUIPMENT COMPANY LTD.

OLD STATION PARK BUILDINGS, ST. JOHN STREET, HORWICH, BOLTON, LANCASHIRE, BL6 7NY, UK